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RPR Test Kit

REF RPR/010* & RPR/012*

*Suffixes indicate change in kit presentation only.

PLASMA TEC

PI-LAB-064/00(01/10/2560)

PRINCIPLE

The RPR (Rapid plasma reagin), test kit is a non-Treponemal test for the qualitative and semi-quantitative detection of syphilis using serum (heated or unheated), or plasma. The RPR test consists of modified VDRL antigen containing carbon particles, which aggregates in the presence of reagin type antibodies in serum or plasma, indicating a positive result.

CLINICAL SIGNIFICANCE

The RPR test is "non-treponemal" in that the antibodies detected are not specific for *T. pallidum*, although their presence in patient's serum or plasma is strongly associated with infection by the organism. This test measures antibody (IgG and IgM) produced in response to lipoidal material released from damaged host cells as well as to lipoprotein-like material released from the spirochaetes. These antibodies tend to disappear after successful treatment of the infection.

WARNINGS AND PRECAUTIONS

For *in vitro* diagnostic use only. For professional use only.

Health and Safety warnings:

All patient samples and reagents should be treated as potentially infectious and the user must wear protective gloves, eye protection and laboratory coats when performing the test.

Non disposable apparatus must be sterilised after use by an appropriate method.

Disposable apparatus must be treated as biohazardous waste and autoclaved or incinerated.

Spillages of potentially infectious material should be absorbed and disposed of as above. The site of spillage must be sterilised with disinfectant or 70% alcohol.

Do not pipette by mouth.

If a reagent vial is compromised, discard the contents immediately.

The product also contains aqueous buffer salts including sodium azide as preservative - see material safety data sheet.

Analytical precautions:

Do not modify the test procedure.

All reagents are ready to use **do not** dilute the reagents in any way.

Reagents and samples to be used at room temperature (18-30°C).

Shake reagents before use to ensure a homogenous suspension.

Avoid cross contamination between different reagents.

Do not interchange droppers. The disposable pipettes should be discarded after a single use.

Discard reagents if they become contaminated or incorrect results are obtained with the controls.

Do not mix reagent of different lots.

STANDARD KIT COMPONENTS

COMPONENT	RPR/010 (100T)	RPR/012 (500T)
RPR Carbon	2ml	10ml
Positive Control	0.5ml	1ml
Negative Control	0.5ml	1ml
Pipette Stirrers	100	500
Dispenser	1	1
Needle	1	1
Reaction Cards	10	50
Instructions For Use	1	1

SPECIMEN AND SAMPLE PREPARATION

This kit can be used with plasma or heated serum. Do not use samples that are haemolysed, turbid or lipaemic. Store at 2-8°C for up to 7 days before testing. If longer storage is required, freeze sample at -20°C or lower. Ensure frozen samples are thawed before testing.

MATERIALS AND EQUIPMENT REQUIRED BUT NOT PROVIDED.

Specimen collection container, Serological pipettes (50 & 100µl), timer, rotating table (100rpm), 0.85% physiological saline, glass test tubes.

RECOMMENDATIONS AND CONTROLS

It is recommended that the positive and negative controls are run with each batch of test specimens. For the assay to be valid the positive control provided should give a strong positive pattern and the negative control provided should give a clearly negative result.

RECOMMENDED PROCEDURE

A. Qualitative Test

1. Draw the sample into the pipstir provided taking care not to transfer and cellular matter.
2. Hold the pipstir vertically above the test card and allow one drop (50µl), of specimen to dispense onto the test card.
3. Spread the specimen evenly over the area of the test circle using the broad end of the pipstir.
4. Ensure the carbon antigen is shaken and homogenous.
5. Attach the needle to the dispenser and withdraw sufficient carbon antigen.
6. Dispense a single drop (20µl) of antigen into the centre of the specimen ensuring the dispenser is held vertically.
7. Rotate card on rotating table for 8 minutes.
8. Read the results macroscopically.

B. Semi Quantative Test

1. Prepare doubling dilutions of the sample from the undiluted specimen to 1:32 using physiological saline.
2. Using the pipstirs, place one drop (50µl), of each dilution onto separate test card circles.
3. Using the broad end of the pipster, spread each dilution over the test circle area from weakest to strongest dilution.
4. Continue procedure as from point 4. of the Qualitative test.

INTERPRETATION OF RESULTS

A. Qualitative Test

Positive: Macroscopic agglutination constitutes a positive test result within the accepted limitations of the test procedure, therefore indicating the presence of the antibodies to *T. Pallidum*.

Negative: No macroscopic agglutination constitutes a negative test result within the accepted limitations of the test procedure, therefore indicating the absence of the antibodies to *T. Pallidum*.

B. Semi Quantative Test

Results can be graded from strong to non-reactive and the titre expressed as the reciprocal of the last dilution showing a positive reaction.

Strong Reactive (SR)	Large clumps of carbon particles with a clear background.
Reactive (R)	Large clumps of carbon particles, more dispersed than strong reactive.
Weak Reactive (WR)	Small clumps of carbon particles with light grey background.
Trace Reactive (TR)	Slight clumping of carbon particles, typically seen as a button of aggregates in the centre of the test circle or dispersed around the edge of the test circle.
Non-Reactive (NR)	A smooth grey pattern or a button of non-aggregated carbon particles in the centre of the test circle.

Reactive samples should be recorded as antibody positive and must be subjected to further tests to determine the presence or absence of specific anti-Treponemal antibody.

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REACTION STABILITY

Read test immediately after rotation. Exercise caution when interpreting results carried out at temperatures other than that recommended.

STORAGE AND SHELF LIFE

Store all reagents upright at 2-8°C.

DO NOT FREEZE THE REAGENT.

Do not use reagents after the stated expiry date.

Discard reagents if they become contaminated.

ALL REAGENTS ARE SUPPLIED READY TO USE

LIMITATIONS OF THE METHOD

1. All reagin tests may give a small proportion of false positive results. Diseases such as infectious mononucleosis, leprosy, lupus erythematosus, vaccinia and viral pneumonia can cause such reactions.
2. Reactive RPR test specimens should be tested with further serological tests (i.e. TPHA and FTA-abs), as with any serological procedure, the diagnosis should not be made on a single reactive result.
3. As with other serological tests, the RPR test cannot distinguish between syphilis and other pathogenic Treponemal infections e.g. Yaws.
4. A diagnosis should always be made in conjunction with clinical findings.
5. False positive/negative results can also occur due to:
 - Contamination of test materials
 - Improper storage, or test temperatures
 - Deviation from the **Recommended Procedures**.

PERFORMANCE CHARACTERISTICS

The following materials were independently tested to compare the performance of the RPR reagent with that of a reference reagent from CDC Atlanta and a reagent from another manufacturer.

Reference plasma panel from CDC Atlanta

Positive test panel from a UK hospital

Normal donor plasma

WHO control

All specimens gave 100% concordant qualitative results with all reagents.

Intra-assay Precision and Accuracy: CV = 0%, accuracy = +/- 0%.

Inter-assay Precision: A standard positive sample tested on eight occasions showed +/- 1 doubling dilution from the nominal value.


DISCLAIMER


The user is responsible for the performance of the reagent by any method other than those mentioned in the **Recommended Procedures**. Any deviations from the **Recommended Procedures** should be validated prior to use.

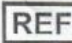BIBLIOGRAPHY

1. McGrew, B.E. et al., *Amer. J. Med Tech.*, **34**, 634 (1968a).
2. McGrew, B.E. et al., *Amer. J. Clin. Path.*, **50**, 52 (1968b).
3. Norins, L.C. *Automation in Clinical Chemistry*, **1**, 157 New York Mediad (1968).
4. Portnoy, J. et al., *U.S. Publ. Health Report*, **77**, 645 (1962).
5. Schroeter, A.L. et al., *Adv. In Automated Analysis*, **1**, 256 N.Y. Mediad
6. Stevens, R.W. and Stroebel, E., *Amer. J. Clin. Path.*, **53**, 32 (1970).
7. Stout, G.W. et al., *J. Conf. Pub. Health Lab. Directors*, **26**, 7, (1968).
8. *Manual Test for Syphilis* PHS Publications No 411, (1969).

TABLE OF SYMBOLS

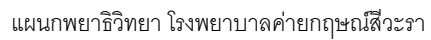SYMBOL	DEFINITION
	Batch Number
	In-vitro Diagnostics
	Catalogue reference
	Store at
	Expiry date
	Manufactured by
	Date of Manufacture
	Read the instructions for use
	Sufficient for

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