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Precision^{14,i,j}

Material	Mean µg/dL [µmol/L]	Standard Deviation (% CV)	
		Repeatability µg/dL [µmol/L]	Within-lab (% CV)
Plasma pool	101 [18.0]	0.5 [0.09] (0.5)	0.7 [0.13] (0.7)
Serum pool 1	95 [16.9]	0.5 [0.09] (0.5)	0.6 [0.11] (0.6)
Serum pool 2	316 [56.6]	1.5 [0.27] (0.5)	3.5 [0.63] (1.1)
Serum pool 3	533 [95.4]	2.4 [0.43] (0.5)	4.2 [0.75] (0.8)
BioRad Lyphochek®			
Control Level 1	231 [41.3]	1.3 [0.23] (0.5)	1.6 [0.29] (0.7)
Control Level 2	50 [8.9]	0.5 [0.09] (1.1)	0.9 [0.16] (1.9)
BioRad Lyphochek®			
Anemia Control Level 1	26 [4.7]	0.3 [0.05] (1.3)	0.5 [0.09] (1.9)
Reduced Sample Volume ^k			
Serum pool 1	103 [18.5]	0.7 [0.13] (0.6)	1.0 [0.18] (0.9)
Serum pool 2	316 [56.6]	1.5 [0.27] (0.5)	3.5 [0.63] (1.1)
Serum pool 3	530 [94.9]	2.9 [0.52] (0.5)	4.2 [0.75] (0.8)
BioRad Lyphochek®			
Anemia Control Level 1	32 [5.7]	0.3 [0.05] (1.3)	0.5 [0.09] (1.9)
BioRad Multiquel®			
Control Level 3	231 [41.3]	1.6 [0.29] (0.7)	2.2 [0.39] (0.9)

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k. Using reduced sample size (25 µL)

Lyphochek® and Multiquel® are registered trademarks of Bio-Rad Laboratories, Irvine, CA 92618.

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Method Comparison

Comparative Method	Regression Statistics ^l		Correlation Coefficient	n ^m
	Slope	Intercept µg/dL [µmol/L]		
Dimension® IRN	0.980	-0.488 [-0.09]	0.9996	147

l. Model equation for regression statistics is: results for Dimension® system = [slope x comparative method results] + intercept.

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n. The range of iron values in the correlation study was 9 to 961 µg/dL [172.0 µmol/L].

Specificity

HIL Interference

The IRON method (using the standard sample size of 40 µL) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-A. Bias is the difference in the results between the control sample (without the interferent) and the test sample (contains the interferent) expressed in percent. Bias exceeding 10% is considered interference.

Substance Tested	Test Concentration SI Units	IRON Concentration µg/mL [µmol/L]	Bias ^o %
Hemoglobin (hemolysate)	50 mg/dL	107 [19.3]	<10
	[0.03 mmol/L] (monomer)		
	200 mg/dL		+10
Bilirubin (unconjugated)	[0.12 mmol/L] (monomer)		
	80 mg/dL	107 [19.3]	<10
	[1368 µmol/L]		
Lipemia (Intralipid®)	3000 mg/dL	128 [22.9]	<10
	[33.9 mmol/L]		

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Substance Tested	Test Concentration SI units	IRON Concentration µg/dL [µmol/L]	Bias ^o %
Hemoglobin (hemolysate)	50 mg/dL	55 [9.8]	<10
	[0.03 mmol/L] (monomer)		
	200 mg/dL	55 [9.8]	+22
Bilirubin (unconjugated)	[0.12 mmol/L] (monomer)		
	80 mg/dL	55 [9.8]	<10
	[1368 µmol/L]		
Lipemia (Intralipid®)	3000 mg/dL	27 [4.8]	<10
	[33.9 mmol/L]		

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Interfering Substances

The IRON method was evaluated for interference according to CLSI/NCCLS EP7-A.¹² Bias is the difference in the results between the control sample (without the interferent) and the test sample (contains the interferent) expressed in mg/dL [$\mu\text{mol/L}$]. Bias exceeding 10% is considered interference.

Interferent	Interferent Concentration SI Units	Iron Concentration $\mu\text{g/dL}$ [$\mu\text{mol/L}$]	Bias $\mu\text{g/dL}$ [$\mu\text{mol/L}$]	Bias (%)
Iron Dextran	60 $\mu\text{g/mL}$	36 [6.4]	+63 [11.3]	+175
	[1074 $\mu\text{mol/L}$]			
Iron Dextran	60 $\mu\text{g/mL}$	131 [23.4]	+69 [12.4]	+53
	[1074 $\mu\text{mol/L}$]			
Hemoglobin	200 mg/dL	55 [9.8]	+12 [2.2]	+22
Hemoglobin	[0.12 mmol/L] (monomer)			
	200 mg/dL	107 [19.2]	+11 [1.9]	+10
Triglycerides (endogenous)	[0.12 mmol/L] (monomer)			
	1109 mg/dL	42 [7.5]	+30 [5.4]	+71 ^{gh}
	[12 mmol/L]			

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Analytical Measurement Range (AMR): 5 – 1000 $\mu\text{g/dL}$ [0.9 – 179.0 $\mu\text{mol/L}$]

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21.1 ใบแทรกน้ำยา IRON Flex® reagent cartridge (PI-LAB-022)

21.2 SOP For Dimension EXL 200 Clinical Chemistry System (MN-LAB-002)

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SIEMENS

REF DF85
PI-LAB-022/00(01/10/2560)

Dimension® clinical chemistry system

Flex® reagent cartridge

IRON

See shaded sections: Updated information from 2008-02 version.

Issue Date 2009-05-15


Iron

Intended Use: The IRON method for the Dimension® clinical chemistry system is an *in vitro* diagnostic test intended to quantitatively measure iron in human serum and plasma. Iron measurements are used in the diagnosis and treatment of diseases such as iron deficiency anemia and other disorders of iron metabolism.¹

Summary: Iron is distributed in the body in such compartments as hemoglobin, tissue, myoglobin, and a labile pool, with the largest amount of iron being found in the hemoglobin of red blood cells or their precursors in bone marrow.¹ Approximately 2.5 mg of physiologic iron is found in plasma, as compared to the approximate 2.5 g of iron contained in hemoglobin.² Disorders of iron metabolism include iron deficiency anemia and iron overload conditions such as hemosiderosis, hemochromatosis, and sideroblastic anemia.³ The IRON method is an adaptation of direct iron assays developed by Smith et al.³ using the chromophore Ferene®. Earlier work by Higgins,⁴ Artiss et al.^{5,6} and Hennessy et al.⁷ demonstrated the high sensitivity of Ferene® and its utility in iron assays. Potential copper interference is minimized by the addition of thiourea.

Ferene® is a registered trademark of Diagnostic Chemicals, LTD., Charlottetown, P.E.I., Canada C1A4H5.

Principles of Procedure: Under acidic conditions, iron (Fe^{3+}) bound to the protein transferrin is released. In the presence of the reducing agent ascorbic acid, (Fe^{3+}) is reduced to (Fe^{2+}). (Fe^{2+}) forms a blue complex with 5,5'-(3-(2-pyridyl)-1,2,4-triazine-5,6-diyl)-bis-2-furansulfonic acid disodium salt (Ferene®). The absorbance of the complex, measured using a bichromatic (600, 700 nm) endpoint technique, is directly proportional to the concentration of iron in the serum.

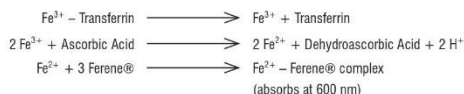Reagents

Wells ^a	Form	Ingredient	Concentration ^b
1 – 4	Liquid	Citric Acid monohydrate	150 mM
		Thiourea Detergent	180 mM
5 – 6 ^c	Liquid	Ferene®	6.0 mM
		Ascorbic Acid	240 mM
		Thiourea	240 mM

a. Wells are numbered consecutively from the wide end of the cartridge.

b. Nominal value per well in a cartridge

c. Wells 5 and 6 contain stabilizer.

Risk and Safety:

Harmful. Contains Thiourea.

R40: Limited evidence of a carcinogenic effect

S36/37: Wear suitable protective clothing and gloves.

Safety data sheets (MSDS/SDS) available on www.siemens.com/diagnostics

Precautions: Used cuvettes contain human body fluids; handle with appropriate care to avoid skin contact or ingestion.

For *in vitro* diagnostic use

Reagent Preparation: All reagents are liquid and ready to use.

Store at: 2 – 8 °C

Expiration: Refer to carton for expiration date of individual unopened reagent cartridges. Sealed wells on the instrument are stable for 30 days.

Open Well Stability: 3 days for wells 1 – 4
14 days for wells 5 – 8

Specimen Collection and Handling: Serum, sodium heparin plasma, and lithium heparin plasma can be collected by normal procedures.⁸ Serum and plasma specimens should be separated from cells within 2 hours after venipuncture.²

Separated specimens are stable for 4 days at room temperature and up to 7 days at 2 – 8 °C. For longer storage, specimens may be frozen at -20 °C for up to 2 months.³ Specimens should be free of particulate matter. To prevent the appearance of fibrin in serum samples, complete clot formation should take place before centrifugation. Clotting time may be increased due to thrombolytic or anticoagulant therapy.

Blood collection tubes containing EDTA, a strong chelator of metal ions, sodium citrate, or a combination of potassium oxalate and sodium fluoride should not be used.⁹

Serum iron concentrations exhibit diurnal variation, with the highest values being obtained in the morning. Because of this diurnal variation, serum iron values may vary by up to 30% during the course of a day.¹⁰ Iron values may remain elevated for several weeks after administration of therapeutic iron-containing compounds such as iron dextran.¹⁰

Hemolyzed samples may give elevated IRON results and should not be used with the IRON method.

Follow the instructions provided with your specimen collection device for use and processing.¹¹

Procedure

Materials Provided
IRON Flex® reagent cartridge, Cat. No. DF85

Materials Required But Not Provided
IRON Calibrator, Cat. No. DC85
Quality Control Materials
Purified Water Diluent (Cat. No. 710615901) or Reagent grade water

Test Steps
Sampling⁴, reagent delivery, mixing, separation, processing and printing of results are automatically performed by the Dimension® system. For details of this processing, refer to your Dimension® Operator's Guide.

d. The sample container must contain sufficient quantity to accommodate the sample volume plus dead volume. Precise container filling is not required.

Test Conditions

	Cuvette 1
Sample Volume	40 µL, 25 µL ^e
Reagent 1 Volume	200 µL
Reagent 2 Volume	70 µL
Temperature	37 °C
Wavelength	600 and 700 nm
Type of Measurement	Bichromatic endpoint

e. An alternate sample size of 25 µL can be programmed; refer to the Dimension® Operator's Guide for the use of an alternate sample size.

Calibration

Assay Range	5 – 1000 µg/dL [0.9 – 179.0 µmol/L] ^f
Calibration Material	Primary standards or secondary calibrators such as IRON Calibrator, Cat. No. DC85
Calibration Scheme	3 levels, n = 3
Units	µg/dL [µmol/L] (µg/dL x 0.179) = [µmol/L]
Typical Calibration Levels	0, 500, 1075 µg/dL [0, 89.5, 192.4 µmol/L]
Calibration Frequency	Every 90 days for any one lot.
A new calibration is required	<ul style="list-style-type: none"> For each new lot of Flex® reagent cartridges After major maintenance or service, if indicated by quality control results As indicated in laboratory quality control procedures When required by government regulations
Assigned Coefficients	Standard sample size = 40 µL C_0 -1.50 C_1 3.46 Alternate sample size = 25 µL C_0 -4.97 C_1 3.49

f. Système International d'Unités [SI Units] are in brackets.

Note: Level 1 calibrator is not included in this carton. Purified Water Diluent (Cat. No. 710615901) or Reagent grade water should be used as the Level 1 calibrator for the IRON method.

Quality Control

At least once each day of use, analyze two levels of a Quality Control (QC) material with known iron concentrations. Follow your laboratory internal QC procedures if the results obtained are outside acceptable limits.

Results: The instrument automatically calculates and prints the concentration of IRON in µg/dL [µmol/L] using the calculation scheme described in your Dimension® Operator's Guide.

Iron determinations can be used in conjunction with Dimension® total iron binding capacity (IBCT) results to calculate percent transferrin saturation (ISAT) and unbound iron binding capacity (UIBC).

Calculated Results: Transferrin Saturation (%): ISAT = 100[IRON/IBCT]
Unbound Iron Binding Capacity: UIBC = [IBCT – IRON]

Refer to your Dimension® Operator's Guide for more information on calculated results.

Results of this test should always be interpreted in conjunction with the patient's medical history, clinical presentation and other findings.

Analytical Measurement Range (AMR): 5 – 1000 µg/dL [0.9 – 179.0 µmol/L]

This is the range of analyte values that can be measured directly from the specimen without any dilution or pretreatment that is not part of the usual analytical process and is equivalent to the assay range.

Samples with results in excess of 1000 µg/dL [179.0 µmol/L] should be repeated on dilution.

Manual Dilution: Dilute with Reagent grade water 1:2 (1 part of sample with 1 part of Reagent grade water) to obtain results within assay range. Enter dilution factor (2). Reassay. Resulting readout is corrected for dilution.

Autodilution (AD): The autodilute sample volume is 20 µL for serum and plasma. Refer to your Dimension® Operator's Guide. Autodilution is not recommended when using the reduced sample size.

Samples with results less than 5 µg/dL [0.9 µmol/L] should be reported as "less than 5 µg/dL [0.9 µmol/L]".

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Limitations of Procedure

The instrument reporting system contains error messages to warn the operator of specific malfunctions. Any report slip containing such error messages should be held for follow-up. Refer to your Dimension® Operator's Guide.

A system malfunction may exist if the following 5-test precision is observed:

Concentration	SD
50 µg/dL [9.0 µmol/L]	>5.1 µg/dL [0.91 µmol/L]
1000 µg/dL [179.0 µmol/L]	>22.1 µg/dL [4.0 µmol/L]

Interfering Substances

The IRON method was evaluated for interference according to CLSI/NCCLS EP7-A.¹² Bias is the difference in the results between the control sample (without the interferent) and the test sample (contains the interferent) expressed in mg/dL [µmol/L]. Bias exceeding 10% is considered interference.

Interferent	Interferent Concentration SI Units	Iron Concentration µg/dL [µmol/L]	Bias µg/dL [µmol/L]	Bias (%)
Iron Dextran	60 µg/mL	36 [6.4]	+63 [11.3]	+175
	[1074 µmol/L]			
Iron Dextran	60 µg/mL	131 [23.4]	+69 [12.4]	+53
	[1074 µmol/L]			
Hemoglobin	200 mg/dL	55 [9.8]	+12 [2.2]	+22
	[0.12 mmol/L] (monomer)			
Hemoglobin	200 mg/dL	107 [19.2]	+11 [1.9]	+10
	[0.12 mmol/L] (monomer)			
Triglycerides (endogenous)	1109 mg/dL	42 [7.5]	+30 [5.4]	+71 ^{a,b}
	[12 mmol/L]			

g. Interference may vary depending upon the lipid composition and subsequent degree of turbidity in the sample.¹³

h. Lipemic samples may be flagged with an "Abnormal Reaction" error message. Refer to your Dimension® Operator's Guide.

Expected Values

Males: 65 – 175 µg/dL [11.6 – 31.3 µmol/L]¹⁰

Females: 50 – 170 µg/dL [9.0 – 30.4 µmol/L]¹⁰

Normal reference intervals can differ by as much as 35% between commercial iron methods;¹ therefore it is advised that each laboratory establish its own expected values for iron as performed on the Dimension® system.

Specific Performance Characteristics

All specific performance characteristic tests were run after normal recommended equipment quality control checks were performed on the Dimension® RxL clinical chemistry system (refer to your Dimension® Operator's Guide).

Precision ^{14,15}				
Material	Mean µg/dL [µmol/L]	Standard Deviation (% CV)		
		Repeatability µg/dL [µmol/L]	Within-lab (% CV)	
Plasma pool	101 [18.0]	0.5 [0.09] (0.5)	0.7 [0.13] (0.7)	
Serum pool 1	95 [16.9]	0.5 [0.09] (0.5)	0.6 [0.11] (0.6)	
Serum pool 2	316 [56.6]	1.5 [0.27] (0.5)	3.5 [0.63] (1.1)	
Serum pool 3	533 [95.4]	2.4 [0.43] (0.5)	4.2 [0.75] (0.8)	
BioRad Lyphocheck®				
Control Level 1	231 [41.3]	1.3 [0.23] (0.5)	1.6 [0.29] (0.7)	
Control Level 2	50 [8.9]	0.5 [0.09] (1.1)	0.9 [0.16] (1.9)	
BioRad Lyphocheck®				
Anemia Control Level 1	26 [4.7]	0.3 [0.05] (1.3)	0.5 [0.09] (1.9)	
Reduced Sample Volume ^x				
Serum pool 1	103 [18.5]	0.7 [0.13] (0.6)	1.0 [0.18] (0.9)	
Serum pool 2	316 [56.6]	1.5 [0.27] (0.5)	3.5 [0.63] (1.1)	
Serum pool 3	530 [94.9]	2.9 [0.52] (0.5)	4.2 [0.75] (0.8)	
BioRad Lyphocheck®				
Anemia Control Level 1	32 [5.7]	0.3 [0.05] (1.3)	0.5 [0.09] (1.9)	
BioRad Multiquel®				
Control Level 3	231 [41.3]	1.6 [0.29] (0.7)	2.2 [0.39] (0.9)	

i. Reproducibility testing was done in accordance with the CLSI/NCCLS Approved Guideline for Evaluation of Precision Performance of Quantitative Measurement Methods (EP5-A2, 2004).

j. Specimens at each level were analyzed in duplicate, twice a day, for 20 days. The repeatability (within-run) and within-lab (total) standard deviations were calculated by analysis of variance method.

k. Using reduced sample size (25 µL)

Lyphocheck® and Multiquel® are registered trademarks of Bio-Rad Laboratories, Irvine, CA 92618.

Method Comparison				
Comparative Method	Regression Statistics ¹			
	Slope	Intercept µg/dL [µmol/L]	Correlation Coefficient	n ^m
Dimension® IRN	0.980	-0.488 [-0.09]	0.9996	147

i. Model equation for regression statistics is: results for Dimension® system = [slope x comparative method results] + intercept.

m. The range of Iron values in the correlation study was: 9 to 963 µg/dL [1.6 to 172.4 µmol/L].

Recommended Sample Types

Matched human serum and plasma samples were analyzed for iron. As shown in the table below, no clinically significant difference was observed when method comparison was conducted using ordinary least squares analysis to fit the regression line between serum and plasma samples.

Sample Type	n ¹	Slope (m)	Intercept (b)	Correlation Coefficient (r)
Sodium heparin plasma versus serum	129	0.988	0.804	0.999
Lithium heparin plasma versus serum	129	0.985	1.42	0.999
Lithium heparin plasma versus sodium heparin plasma	129	0.997	0.666	0.999

n. The range of Iron values in the correlation study was 9 to 961 µg/dL [172.0 µmol/L].

Specificity

HIL Interference
The IRON method (using the standard sample size of 40 µL) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-A. Bias is the difference in the results between the control sample (without the interferent) and the test sample (contains the interferent) expressed in percent. Bias exceeding 10% is considered interference.

Substance Tested	Test Concentration SI Units	IRON Concentration µg/dL [µmol/L]	Bias ^a %
Hemoglobin (hemolysate)	50 mg/dL	107 [19.3]	<10
	[0.03 mmol/L] (monomer)		
Bilirubin (unconjugated)	200 mg/dL		+10
	[0.12 mmol/L] (monomer)		
Lipemia (Intralipid®)	80 mg/dL	107 [19.3]	<10
	[1368 µmol/L]		
Lipemia (Intralipid®)	3000 mg/dL	128 [22.9]	<10
	[33.9 mmol/L]		

Intralipid® is a registered trademark of Fresenius Kabi AG, Bad Homburg, Germany.

The IRON method (using the reduced sample size of 25 µL) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-A. Bias is the difference in the results between the control sample (without the interferent) and the test sample (contains the interferent) expressed in percent. Bias exceeding 10% is considered interference.

Substance Tested	Test Concentration SI Units	IRON Concentration µg/dL [µmol/L]	Bias ^a %
Hemoglobin (hemolysate)	50 mg/dL	55 [9.8]	<10
	[0.03 mmol/L] (monomer)		
Bilirubin (unconjugated)	200 mg/dL	55 [9.8]	+22
	[0.12 mmol/L] (monomer)		
Lipemia (Intralipid®)	80 mg/dL	55 [9.8]	<10
	[1368 µmol/L]		
Lipemia (Intralipid®)	3000 mg/dL	27 [4.8]	<10
	[33.9 mmol/L]		

o. Analyte results should not be corrected based on this bias.

Non-Interfering Substances

The following substances do not interfere with the IRON method when present in serum in the amounts indicated. Systematic inaccuracies (bias) due to these substances are less than 10% at an iron concentration of 26 to 38 µg/dL [4.7 to 6.8 µmol/L] and 118 to 136 µg/dL [21.1 to 24.3 µmol/L].

Substance	Test Concentration	SI Units
Acarbose	180 µg/mL	200 µmol/L
Acetaminophen	20 mg/dL	1323 µmol/L
Allopurinol	40 µg/dL	253 µmol/L
Amikacin	15 mg/dL	256 µmol/L
Amiodorone HCl	6 µg/mL	8.8 µmol/L
Ampicillin	5.3 mg/dL	152 µmol/L
Ascorbic Acid	5 mg/dL	227 µmol/L
Atenolol	1 mg/dL	37.5 µmol/L
Atorvastatin	600 µg/L	0.4 µmol/L
Caffeine	10 mg/dL	515 µmol/L
Calcitriol	0.3 µg/mL	0.7 µmol/L
Calcium	15 mg/dL	1.4 mmol/L
Captopril	22 µg/mL	101 µmol/L
Carbamazepine	12 mg/dL	508 µmol/L
Chloramphenicol	25 mg/dL	774 µmol/L
Chlordiazepoxide	2 mg/dL	66 µmol/L
Chlorpromazine	5 mg/dL	157 µmol/L
Cholesterol	500 mg/dL	12.9 mmol/L
Cinacalcet hydrochloride	18 ng/mL	46 nmol/L
Cimetidine	10 mg/dL	396 µmol/L
Copper	300 µg/dL	22.3 µmol/L
Creatinine	30 mg/dL	2652 µmol/L
Deferoxamine	250 ng/dL	3.8 nmol/L
Dextran 40	6000 mg/dL	1500 µmol/L
Diazepam	2 mg/dL	70 µmol/L
Diliazem hydrochloride	40 ng/mL	89 nmol/L
Digoxin	5 ng/mL	6.2 nmol/L
Disopyramide phosphate	4 mg/dL	91.8 µmol/L
Epoetin alfa	456 mU/mL	456 U/L
Erythromycin	20 mg/dL	273 µmol/L
Ethanol	400 mg/dL	87 mmol/L
Ethosuximide	30 mg/dL	2125 µmol/L
Fenofibrate	45 µg/dL	125 µmol/L
Ferritin	200 ng/mL	449 pmol/L
Folic acid	1 ng/mL	2.3 nmol/L
Fluvastatin	48 mg/dL	1.1 mmol/L
Furosemide	6 mg/dL	181 µmol/L
Gemfibrozil	75 µg/mL	299 µmol/L
Gentamicin	12 mg/dL	251 µmol/L
Glyburide	2 µg/mL	4 µmol/L
Heparin	8 U/mL	8000 U/L
Hydrochlorothiazide	5.9 µg/mL	198 µmol/L
Ibuprofen	50 mg/dL	2425 µmol/L
Immunoglobulin G (IgG)	5 g/dL	50 g/L
Insulin	0.018 U/mL	18 U/L
Isosorbide dinitrate	150 ng/mL	635 µmol/L
Lidocaine	6 mg/dL	256 µmol/L
Lithium	3.5 mg/dL	5043 µmol/L
Losartan potassium	10 mg/dL	0.22 mmol/L
Magnesium	15 mg/dL	1.6 mmol/L
Metformin	40 µg/mL	241 µmol/L
Nateglinide	72 µg/mL	227 µmol/L
Niacin	1.2 mg/mL	9.7 mmol/L
Nicotine	2 mg/dL	123 µmol/L
Nitrofurantoin	4 µg/mL	17 µmol/L
Nortryptiline	1 µg/mL	3 µmol/L

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Analytical Sensitivity: 5 µg/dL [0.9 µmol/L]


The analytical sensitivity represents the lowest concentration of iron that can be distinguished from zero. This sensitivity is defined as the mean value ($n = 20$) plus two standard deviations of the low level (0 µg/dL [0 µmol/L]) IRON Calibrator, Purified Water Diluent (Cat. No. 710615901) or Reagent grade water should be used as the Level 1 calibrator.

Symbols Key: See adjacent panel.

Bibliography: See adjacent panel.

Dimension® and Flex® are trademarks of Siemens Healthcare Diagnostics Inc.

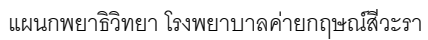©2009 Siemens Healthcare Diagnostics Inc.
All rights reserved.

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