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i. Model equation for regression statistics is: Results of Dimension® TP system
= [slope x (comparative method results)] + intercept.

j. Range of Samples: 2.8 – 8.8 g/dL [28 – 88 g/L].

k. Range of Samples: 3.0 – 11.6 g/dL [30 – 116 g/L].

l. Model equation for regression statistics is: Results of Dimension® TP system using reduced sample size
(10 µL) = [slope x Results of Dimension® TP system using standard sample size (15 µL)] + intercept.

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Specificity

HIL Interference			
The TP method (using the standard sample size of 15 µL) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-P. Bias, defined as the difference between the control sample (does not contain interferent) and the test sample (contains the interferent), is shown in the table below. Bias exceeding 10% is considered "interference".			
Substance Tested	Test Concentration SI Units	TP Concentration g/dL [g/L]	Bias* %
Hemoglobin (hemolysate)	300 mg/dL [0.19 mmol/L] (monomer)	3.9 [39]	<10
Bilirubin (unconjugated)	5 mg/dL [86 µmol/L]	3.8 [38]	<10
Lipemia (Intralipid®)	200 mg/dL [2.26 mmol/L]	3.9 [39]	<10
	600 mg/dL [6.78 mmol/L]	3.9 [39]	n
The TP method (using the reduced sample size of 10 µL) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-P. Bias, defined as the difference between the control sample (does not contain interferent) and the test sample (contains the interferent), is shown in the table below. Bias exceeding 10% is considered "interference".			
Substance Tested	Test Concentration SI units	TP Concentration g/dL [g/L]	Bias %
Hemoglobin (hemolysate)	200 mg/dL [0.12 mmol/L] (monomer)	3.2 [32]	<10
Bilirubin (unconjugated)	5 mg/dL [86 µmol/L]	3.2 [32]	<10
Lipemia (Intralipid®)	200 mg/dL [2.26 mmol/L]	3.1 [31]	<10
	600 mg/dL [6.78 mmol/L]	3.2 [32]	n
n. The interference testing at this level tripped a test report message; therefore the magnitude of the interference could not be determined.			
* Analyte results should not be corrected based on this bias.			

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21.2 SOP For Dimension EXL 200 Clinical Chemistry System (MN-LAB-002)

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REF DF73

PI-LAB-014/00(01/10/2560)

Dimension® clinical chemistry system

Flex® reagent cartridge

TP

See shaded sections: Updated information from 2003-04 version.

Issue Date 2008-04-23

Total Protein


Intended Use: The TP method used on the Dimension® clinical chemistry system is an *in vitro* diagnostic test intended for the quantitative determination of total protein in human serum and heparinized plasma. Measurements of total protein are used in the diagnosis and treatment of a variety of diseases involving the liver, kidney or bone marrow as well as metabolic or nutritional disorders.

Summary: The total protein method is a modification of the biuret reaction first introduced by Kingsley¹ and later modified by Henry² and presented as the method of choice for serum by Henry.³ This method incorporates tartrate as a complexing agent to prevent precipitation of Cu(OH)₂. Serum blanking increases method sensitivity and minimizes spectral interference from lipemia.

Principles of Procedure: Cupric ion (Cu²⁺) reacts with the peptide linkages (-C-NH-CH-C-NH-) of protein in a basic solution.

The blue copper (II) protein complex thus formed is proportional to the total protein concentration in the sample and is measured using a bichromatic (540, 700 nm) endpoint technique.

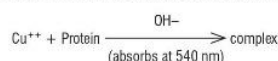Reagents

Wells ^a	Form	Ingredient	Concentration ^b
1 – 3	Liquid	Potassium Sodium Tartrate NaOH	1.089 g/mL
4 – 6	Liquid	Cupric Sulfate	0.015 mol/L

a. Wells are numbered consecutively from the wide end of the cartridge.
b. Nominal value per test at manufacture.

Risk and Safety:

R52/53: Harmful to aquatic organisms; may cause long-term adverse effects in the aquatic environment.

Safety data sheets (MSDS/SDS) available on www.siemens.com/diagnostics

Precautions: Used cuvettes contain human body fluids; handle with appropriate care to avoid skin contact and ingestion.

For *in vitro* diagnostic use

Reagent Preparation: All reagents are liquid and ready to use.

Store at: 2 – 8 °C

Expiration: Refer to shelf carton for expiration date of individual unopened reagent cartridges. Sealed cartridge wells on the instrument are stable for 30 days.

Open Well Stability: 5 days for wells 1 – 6

Specimen Collection and Handling: Serum and plasma can be collected using recommended procedures for collection of diagnostic blood specimens by venipuncture.^{4,5}

Follow the instructions provided with your specimen collection device for use and processing.⁶

Complete clot formation should take place before centrifugation.⁷

Specimens should be free of particulate matter.

Separated specimens are stable for 8 hours at room temperature, 72 hours at 2 – 8 °C and 6 months when frozen at -20 °C or colder.^{1,8}

Corvac® and SST® collection tubes do not affect the TP method.

Blood collection tubes containing EDTA, lithium heparin, potassium oxalate or sodium fluoride do not affect the TP method.

Corvac® is a registered trademark of Monoject, Division of Sherwood Medical, St. Louis, MO.

SST® is a registered trademark of Becton-Dickinson, Rutherford, NJ.

Procedure

Materials Provided	
TP Flex® reagent cartridge, Cat. No. DF73	
Materials Required But Not Provided	
Total Protein/Albumin Calibrator, Cat. No. DC31	
Quality Control Materials	

Test Steps

Sampling,⁹ reagent delivery, mixing, processing, and printing of results are automatically performed by the Dimension® system. For details of this processing, refer to your Dimension® Operator's Guide.

c. The sample container (if not a primary tube) must contain sufficient quantity to accommodate the sample volume plus dead volume. Precise container filling is not required.

Test Conditions

Sample Size	15 µL, (10 µL) ^d
Reagent 1 Volume	85 µL
Reagent 2 Volume	85 µL
Diluent Volume	315 µL
Temperature	37 °C
Wavelength	540 and 700 nm
Type of Measurement	Bichromatic endpoint

d. An alternate sample size (10 µL) can be programmed; refer to the Operator's Guide for the use of alternate sample size.

Calibration

Assay Range	2.0 – 12.0 g/dL [20 – 120 g/L]*
Calibration Material	Total Protein/Albumin Calibrator, Cat. No. DC31
Calibration Scheme	3 levels, n = 3
Units	g/dL [g/L] (g/dL x 10) = [g/L]
Typical Calibration Levels	2.0, 6.0, 10.0 g/dL [20, 60, 100 g/L]
Calibration Frequency	Every 3 months for any one lot
A new calibration is required	<ul style="list-style-type: none"> For each new lot of Flex® reagent cartridges After major maintenance or service, if indicated by quality control results As indicated in laboratory quality control procedures When required by government regulations
Assigned Coefficients	Standard sample size = 15 µL C ₁ 3.700 C ₂ 0.0222 Alternate sample size = 10 µL C ₁ 5.357 C ₂ 0.031

e. Système International d'Unités [SI Units] are in brackets.

Quality Control

At least once each day of use, analyze two levels of a Quality Control (QC) material with known total protein concentrations.

Follow your laboratory internal QC procedures if the results obtained are outside acceptable limits.

Results: The instrument automatically calculates and prints the concentration of total protein in g/dL [g/L] using the calculation scheme illustrated in your Dimension® Operator's Guide.

Results of this test should always be interpreted in conjunction with the patient's medical history, clinical presentation and other findings.

Analytical Measurement Range (AMR): 2.0 – 12.0 g/dL [20 – 120 g/L]

This is the range of analyte values that can be directly measured on the specimen without any dilution or pretreatment that is not part of the usual analytical process and is equivalent to the assay range.

Samples with results in excess of 12.0 g/dL [120 g/L] should be repeated on dilution.

Manual Dilution: Make appropriate dilution with Reagent grade water to obtain result within the assay range. Enter dilution factor.

Reassay. Resulting readout is corrected for dilution.

Autodilution (AD): Refer to your Dimension® Operator's Guide.

Results less than 2.0 g/dL [20 g/L] should be reported as "less than 2.0 g/dL [20 g/L]" instead of the numerical value.

Limitations of Procedure

The instrument reporting system contains error messages to warn the operator of specific malfunctions. Any report slip containing such error messages should be held for follow-up. Refer to your Dimension® Operator's Guide.

A system malfunction may exist if the following 5-test precision is observed at the standard sample size (15 µL):

Concentration	SD
6.8 g/dL [68 g/L]	>0.12 g/dL [1.2 g/L]

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Interfering Substances

Dextran 40 of 1500 mg/dL [375 µmol/L] increases a TP result of 7.0 g/dL [70 g/L] by 17%.
Immunoglobulin G of 2.5 g/dL [25 g/L] increases a TP result of 7.0 g/dL [70 g/L] by 25%.

Using standard sample size (15 µL):

Hemoglobin (hemolysate) of 500 mg/dL [0.31 mmol/L] (monomer) increases a TP result of 3.9 g/dL [39 g/L] by 11%.

Bilirubin (unconjugated) of 20 mg/dL [342 µmol/L] decreases a TP result of 3.8 mg/dL [38 g/L] by -11%.

Lipemia (Intralipid®) of 600 mg/dL [6.78 mmol/L] and above tripped an error flag on this method, so the magnitude of the interference is not available.

For reduced sample size (10 µL):

Hemoglobin (hemolysate) of 300 mg/dL [0.19 mmol/L] (monomer) increases a TP result of 3.2 g/dL [32 g/L] by 11%.

Bilirubin (unconjugated) of 20 mg/dL [342 µmol/L] decreases a TP result of 3.2 mg/dL [32 g/L] by -13%.

Lipemia (Intralipid®) of 600 mg/dL [6.78 mmol/L] and above tripped an error flag on this method, so the magnitude of the interference is not available.

Intralipid® is a registered trademark of Fresenius Kabi AG, Bad Homburg, Germany.

Expected Values^h: 6.4 – 8.2 g/dL [64 – 82 g/L]

This reference population consisted of the following:

118 males, ages 20 – 65

108 females, ages 20 – 65

The reference interval was calculated non-parametrically and represents the central 95% of the population, using serum as the sample.

A reference interval for plasma would be 0.3 g/dL higher due to the presence of fibrinogen in the plasma.

Each laboratory should establish its own reference interval for total protein as performed on the Dimension® system.

Specific Performance Characteristics^f

Material	Mean g/dL [g/L]	Precision ^g	
		Standard Deviation (% CV) Within-run	Total
Multiquel® Chemistry Control			
Normal	6.06 [61]	0.04 [0.4] (0.7)	0.09 [0.9] (1.5)
Abnormal	4.64 [46]	0.05 [0.5] (1.1)	0.10 [1.0] (2.2)
Moni-Trol® ^h			
Level 1	7.34 [73]	0.04 [0.4] (0.6)	0.06 [0.6] (0.8)
Level 2	4.80 [48]	0.05 [0.5] (0.9)	0.07 [0.7] (1.4)

f. All specific performance characteristics tests were run after normal recommended equipment quality control checks were performed (refer to your Dimension® Operator's Guide).

g. Specimens at each level were analyzed in triplicate for 20 runs. The within-run and total standard deviations were calculated by the analysis of variance method.

h. Using reduced sample size (10 µL).

Multiquel® Chemistry Control is a registered trademark of Bio-Rad Laboratories, Inc. Irvine, CA 92618.

Moni-Trol® is a registered trademark of Medical Analysis Systems Inc., Camarillo, CA 93012-8058.

Method Comparison

Regression Statistics

Comparative Method	Slope	Intercept g/dL [g/L]	Correlation Coefficient	n
TP Method on the aca® discrete clinical analyzer ⁱ	0.984	0.04 [0.4]	0.966	60 ^j
Serum TP method on the Dimension® vs. plasma	0.982	0.37 [4.9]	0.9877	52 ^k
Reduced vs. Standard sample size ^l	0.975	0.12 [1.2]	0.996	60 ^m

i. Model equation for regression statistics is: Results of Dimension® TP system = [slope x (comparative method results)] + intercept.

j. Range of Samples: 2.8 – 8.8 g/dL [28 – 88 g/L].

k. Range of Samples: 3.0 – 11.6 g/dL [30 – 116 g/L].

l. Model equation for regression statistics is: Results of Dimension® TP system using reduced sample size (10 µL) = [slope x Results of Dimension® TP system using standard sample size (15 µL)] + intercept.

m. Range of samples: 3.9 – 9.7 g/dL [39 – 97 g/L].

Specificity

HIL Interference

The TP method (using the standard sample size of 15 µL) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-P Bias, defined as the difference between the control sample (does not contain interferent) and the test sample (contains the interferent), is shown in the table below. Bias exceeding 10% is considered "interference".

Substance Tested	Test Concentration SI Units	TP Concentration g/dL [g/L]	Bias* %
Hemoglobin (hemolysate)	300 mg/dL [0.19 mmol/L] (monomer)	3.9 [39]	<10
Bilirubin (unconjugated)	5 mg/dL [86 µmol/L]	3.8 [38]	<10
Lipemia (Intralipid®)	200 mg/dL [2.26 mmol/L]	3.9 [39]	<10
	600 mg/dL [6.78 mmol/L]	3.9 [39]	n

The TP method (using the reduced sample size of 10 µL) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-P Bias, defined as the difference between the control sample (does not contain interferent) and the test sample (contains the interferent), is shown in the table below. Bias exceeding 10% is considered "interference".

Substance Tested	Test Concentration SI units	TP Concentration g/dL [g/L]	Bias %
Hemoglobin (hemolysate)	200 mg/dL [0.12 mmol/L] (monomer)	3.2 [32]	<10
Bilirubin (unconjugated)	5 mg/dL [86 µmol/L]	3.2 [32]	<10
Lipemia (Intralipid®)	200 mg/dL [2.26 mmol/L]	3.1 [31]	<10
	600 mg/dL [6.78 mmol/L]	3.2 [32]	n

n. The interference tested at this level tripped a test report message; therefore the magnitude of the interference could not be determined.

* Analyte results should not be corrected based on this bias.

Non-Interfering Substances

The following substances have a negligible effect (< 0.2 g/dL [< 2 g/L]) on the TP method at the concentrations indicated:

Substance	Test Concentration	SI Units
Acetaminophen	0.025 mg/dL	1.66 µmol/L
Amikacin	15 mg/dL	256 µmol/L
Ampicillin	5.3 mg/dL	152 µmol/L
Ascorbic Acid	5 mg/dL	227 µmol/L
Caffeine	6 mg/dL	308 µmol/L
Carbamazepine	3 mg/dL	127 µmol/L
Chloramphenicol	5 mg/dL	155 µmol/L
Chlordiazepoxide	1 mg/dL	33.3 µmol/L
Chlorpromazine	0.2 mg/dL	6.27 µmol/L
Cholesterol	500 mg/dL	12.9 mmol/L
Cimetidine	2 mg/dL	79.2 µmol/L
Creatinine	30 mg/dL	2652 µmol/L
Diazepam	0.5 mg/dL	17.6 µmol/L
Erythromycin	6 mg/dL	81.6 µmol/L
Ethanol	400 mg/dL	86.8 mmol/L
Ethosuximide	25 mg/dL	1770 µmol/L
Furosemide	6 mg/dL	181 µmol/L
Gentamicin	12 mg/dL	251 µmol/L
Heparin	3 U/mL	3000 U/L
Ibuprofen	50 mg/dL	2425 µmol/L
Lidocaine	1.2 mg/dL	51.2 µmol/L
Lithium	2.3 mg/dL	3.2 mmol/L
Nicotine	0.1 mg/dL	6.2 µmol/L
Penicillin G	25 U/mL	25000 U/L
Pentobarbital	8 mg/dL	354 µmol/L
Phenobarbital	10 mg/dL	421 µmol/L
Phenytoin	5 mg/dL	198 µmol/L
Primidone	4 mg/dL	183 µmol/L
Propoxyphene	0.2 mg/dL	4.91 µmol/L
Salicylic Acid	60 mg/dL	4.34 mmol/L
Theophylline	4 mg/dL	222 µmol/L
Urea	500 mg/dL	83.3 mmol/L
Uric Acid	20 mg/dL	1190 µmol/L
Valproic Acid	50 mg/dL	3467 µmol/L


Analytical Sensitivity: ≤ 2.0 mg/dL [20 g/L]

The analytical sensitivity represents the lowest concentration of TP that can be distinguished from zero. This sensitivity is defined as the low assay range of TP.

Symbols Key: See adjacent panel.

Dimension®, aca® and Flex® are trademarks of Siemens Healthcare Diagnostics Inc.

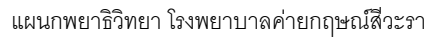©2008 Siemens Healthcare Diagnostics Inc.
All rights reserved.

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