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
- 1.1 Bone diseases
- 1.2 Chronic renal diseases ,dialysis patients
- 1.3 Status post thyroid surgery
- 1.4 Parathyroid disease
- 1.5 Patients with nephrolithiasis
- 1.6 Chronic alcoholism
- 1.7 Patients requiring intensive medical care (parenteral nutrition, mechanical ventilation)
- 1.8 Suspected vitamin D deficiency (malabsorption syndrome)
- 1.9 Muscle weakness, bone pain

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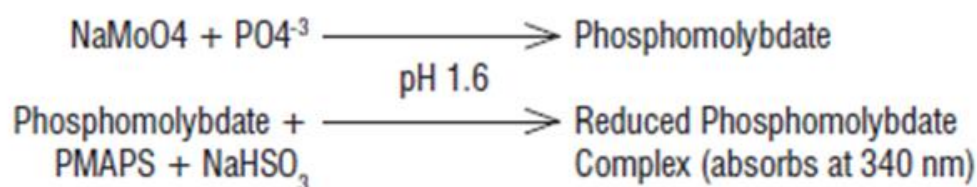Specific Performance Characteristics*

Material	Precision ^f		
	Mean mg/dL [mmol/L]	Standard Deviation (% CV)	
		Within-run	Between-day
QCS™ Control Serum Assayed			
Normal	3.6 [1.16]	0.05 [0.02] (1.4)	0.13 [0.04] (3.6)
Abnormal	8.5 [2.74]	0.08 [0.03] (0.9)	0.17 [0.05] (2.0)
Fisher Urine Chemistry Control			
Level 1	61.5 [19.8]	1.1 [0.3] (1.8)	7.3 [2.3] (11.9)
Level 2	70.3 [22.7]	1.3 [0.4] (1.8)	7.3 [2.3] (10.4)
Urine Pool	43.9 [14.2]	0.6 [0.2] (1.4)	3.6 [1.2] (8.2)

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QCS™ is a trademark of Gilford, Irvine, CA 92714.

Method Comparison

Regression Statistics^g

Comparative Method	Slope	Intercept mg/dL [mmol/L]	Correlation Coefficient	n
PHOS method on the				
aca® discrete clinical analyzer				
Serum	1.00	0.0	0.976	112
Urine	1.00	-0.05 [-0.02]	0.992	122

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Substance Tested	Test Concentration SI Units	PHOS Concentration mg/dL [mmol/L]	Bias % ^h
Hemoglobin (hemolysate)	50 mg/dL [0.03 mmol/L] (monomer)	1.1 [0.36 L]	<10
Bilirubin (unconjugated)	5 mg/dL [86 µmol/L]	1.0 [0.32]	<10
Lipemia (Intralipid®)	200 mg/dL [2.29 mmol/L]	1.1 [0.36]	<10

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21.1 ใบแทรกน้ำยา PHOS Flex® reagent cartridge (PI-LAB-012)

21.2 SOP For Dimension EXL 200 Clinical Chemistry System (MN-LAB-002)

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Flex® reagent cartridge

PHOS

See shaded sections: Updated information from 2001-01 version.


Issue Date 2008-03-27

Phosphorus

Intended Use: The PHOS method used on the Dimension® clinical chemistry system is an *in vitro* diagnostic test intended for the quantitative determination of phosphorus in human serum, plasma and urine.

Summary: The phosphorus (PHOS) method is a modification of the classical phosphomolybdate method introduced by Fiske and Subbarow¹ and uses a mixture of p-methylaminophenol sulfate and bisulfite to reduce the phosphomolybdate as reported by Gomori² and by Drewes.³ The PHOS method measures the absorption of the reduced phosphomolybdate complex.⁴ Interference caused by protein precipitation is eliminated using the solubilizing agent, lithium dodecyl sulfate. Because of the pre-blanking of the sample, bilirubin interference is considerably minimized.

Principles of Procedure: Inorganic phosphate combines with molybdate (MoO_4) in an acid solution to form a complex which is reduced by p-methylaminophenol sulfate (PMAPS) and bisulfite. The 340 nm absorbance of the reduced phosphomolybdate solution is proportional to the inorganic phosphorus concentration and is measured using a bichromatic endpoint technique.

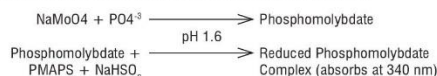Reagents

Wells ^a	Form	Ingredient	Concentration ^b
1 – 6	Liquid	PMAPS	0.64 mmol/L
		Sodium bisulfite	8.55 mmol/L
7	Liquid	Lithium Dodecyl Sulfate	6.15 mmol/L
8	Liquid	Sodium Molybdate	1.10 mmol/L
		H_2SO_4	

a. Wells are numbered consecutively from the wide end of the cartridge.
b. Nominal value per test at manufacture.

Risk and Safety:

Irritant. Contains mixture of 5-chloro-2-methyl-2H-isothiazol-3-one and 2-methyl-2H-isothiazol-3-one (3:1).
R36/38: Irritating to eyes and skin.
R52/53: Harmful to aquatic organisms, may cause long-term adverse effects in the aquatic environment.
S37/39: Wear suitable gloves and eye/face protection.

Safety data sheets (MSDS/SDS) available on www.siemens.com/diagnostics

Precautions: Used cuvettes contain human body fluids; handle with appropriate care to avoid skin contact and ingestion.

For *in vitro* diagnostic use

Reagent Preparation: All reagents are liquid and ready to use.

Store at: 2 – 8 °C

Expiration: Refer to carton for expiration date of individual unopened reagent cartridges. Sealed cartridge wells on the instrument are stable for 30 days.

Open Well Stability: 3 days for wells 1 – 6
30 days for wells 7 – 8

Specimen Collection and Handling: Serum and plasma can be collected using recommended procedures for collection of diagnostic blood specimens by venipuncture.⁵

Follow the instructions provided with your specimen collection device for use and processing.⁶

Complete clot formation should take place before centrifugation. Serum or plasma should be separated from red cells within 1 hour because erythrocytes contain phosphate concentrations several times greater than those found in the serum.^{7,8}

Specimens are stable for 8 hours at room temperature, 2 days at 2 – 8 °C. For longer storage, specimens may be frozen at -20 °C or colder.⁷

Corvac® and SST® collection tubes do not affect the PHOS method.

Collection of 24 hour urine specimen should be made in a container with 10 to 20 mL of 6M HCl. If preservatives are not used during 24 hour urine collection, acidify the sample to below pH 3.0 before analysis.⁹

Corvac® is a registered trademark of Monoject, Division of Sherwood Medical, St. Louis, MO.
SST® is a registered trademark of Becton-Dickinson, Rutherford, NJ.

Procedure

Materials Provided

PHOS Flex® reagent cartridge, Cat. No. DF61

Materials Required But Not Provided

CHEM II Calibrator, Cat. No. DC20
Quality Control Materials

Test Steps

Sampling,² reagent delivery, mixing, processing and printing of results are automatically performed by the Dimension® system. For details of this processing, refer to your Dimension® Operator's Guide.

c. The sample container (if not a primary tube) must contain sufficient quantity to accommodate the sample volume plus the dead volume; precise container filling is not required.

Test Conditions

Sample Size	3 µL
Reagent 1 Volume	50 µL
Reagent 2 Volume	20 µL
Reagent 3 Volume	20 µL
Diluent Volume	350 µL
Temperature	37 °C
Wavelength	*340 – 383 nm 340 – 700 nm
Type of Measurement	Bichromatic endpoint

*Dimension® AR 4.4 uses the 340 – 383 nm wavelength pair.

Calibration

Assay Range	0 – 9.0 mg/dL [0 – 2.90 mmol/L]
Calibration Material	CHEM II Calibrator, Cat. No. DC20
Calibration Scheme	3 Levels, n = 3
Units	mg/dL [mmol/L] ^a (mg/dL x 0.323) = [mmol/L]
Typical Calibration Levels	2.0, 5.0, 8.0 mg/dL [0.65, 1.62, 2.58 mmol/L]
Calibration Frequency	Every 3 months for any one lot.
A new calibration is required	<ul style="list-style-type: none"> For each new lot of Flex® reagent cartridges After major maintenance or service, if indicated by quality control results As indicated in laboratory quality control procedures When required by government regulations
Assigned Coefficients	C ₀ -1.100 C ₁ 0.077

d. Système International d'Unités [SI Units] are in brackets.

Quality Control

At least once each day of use, analyze two levels of a Quality Control (QC) material with known phosphorus concentrations. Follow your laboratory internal QC procedures if the results obtained are outside acceptable limits.

Results: The instrument automatically calculates and prints the concentration of PHOS in mg/dL [mmol/L] using the calculation scheme illustrated in your Dimension® Operator's Guide.

Results of this test should always be interpreted in conjunction with the patient's medical history, clinical presentation and other findings.

Analytical Measurement Range (AMR): 0 – 9.0 mg/dL [0 – 2.90 mmol/L]

This is the range of analyte values that can be directly measured on the specimen without any dilution or pretreatment that is not part of the usual analytical process and is equivalent to the assay range.

Samples with results in excess of 9.0 mg/dL [2.90 mmol/L] should be repeated on dilution.

Manual Dilution: Serum/Plasma: Dilute with Reagent grade water to obtain result within the assay range. Enter dilution factor.
Urine: Dilute 1 part urine: 9 parts Reagent grade water or equivalent. Enter dilution factor of 10.
Serum/Plasma/Urine: Reassay. Resulting readout is corrected for dilution.

Autodilution (AD):

(for serum, plasma): Refer to your Dimension® Operator's Guide.

Automated Urine Dilution:

(AUD) (for urine): Refer to your Dimension® Operator's Guide.

Limitations of Procedure

The instrument reporting system contains error messages to warn the operator of specific malfunctions. Any report slip containing such error messages should be held for follow-up. Refer to your Dimension® Operator's Guide.

A system malfunction may exist if the following 5-test precision is observed:

Concentration	SD
2.0 mg/dL [0.65 mmol/L]	>0.15 mg/dL [0.05 mmol/L]
8.0 mg/dL [2.58 mmol/L]	>0.40 mg/dL [0.13 mmol/L]

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Interfering Substances

Hemolyzed samples may give spuriously elevated phosphorus results. Bias from hemolysis may result from inorganic phosphates produced by the action of phosphatases on organic phosphates, both of which are released from red blood cells upon hemolysis.¹⁰

Hemoglobin (hemolysate) at 200 mg/dL [0.12 mmol/L] (monomer) increases a PHOS result of 1.1 mg/dL [0.36 mmol/L] by 27%.

Bilirubin at 20 mg/dL [342 mmol/L] decreases a PHOS result of 1.0 mg/dL [0.4 mmol/L] by 21%.

Lipemia (Intralipid®) at 600 mg/dL [6.8 mmol/L] decreases a PHOS result of 1.0 mg/dL [0.4 mmol/L] by 16%.

Serum or plasma should be separated from erythrocytes soon after the specimen is collected; otherwise, the phosphatase and inorganic phosphates present in the red cells will interfere and give falsely high values.

Creatine phosphate and phosphoenol pyruvate liberate phosphorus under the reaction conditions of this method and can thereby increase results.

Mannitol at concentrations of 500 mg/dL [27.4 mmol/L] or greater causes a decrease of 10% or more.

Lipemia with triglyceride of 600 mg/dL [6.86 mmol/L] and bilirubin of 20 mg/dL [342 µmol/L] increase the PHOS result by 0.2 mg/dL [0.06 mmol/L] and 0.8 mg/dL [0.26 mmol/L] respectively, at a phosphorus concentration of 6.9 mg/dL [2.23 mmol/L].

For Dimension® systems using 340 – 383 nm, hemoglobin of 500 mg/dL [0.31 mmol/L] decreases the PHOS result by 1.2 mg/dL [0.39 mmol/L] at a PHOS concentration of 6.5 mg/dL [2.10 mmol/L].

For Dimension® systems using 340 – 700 nm, interference from hemoglobin at 1000 mg/dL [0.62 mmol/L] was less than 10% at a PHOS concentration of 4.1 mg/dL [1.33 mmol/L].

Intralipid® is a registered trademark of Fresenius Kabi AG, Bad Homburg, Germany.

Expected Values

Serum: 2.5 – 4.9 mg/dL [0.81 – 1.58 mmol/L]¹¹

Urine: 0.4 – 1.3 g/24 hr [12 – 42 mmol/24 hr]⁹

This reference population consisted of 136 males and 130 females, ages 17 – 60.

Each reference interval was calculated non-parametrically and represents the central 95% of the population.

Each laboratory should establish its own reference intervals for PHOS as performed on the Dimension® system.

Specific Performance Characteristics^a

Material	Precision ^a		
	Mean mg/dL [mmol/L]	Standard Deviation (% CV)	
		Within-run	Between-day
QCS™ Control Serum Assayed			
Normal	3.6 [1.16]	0.05 [0.02] (1.4)	0.13 [0.04] (3.6)
Abnormal	8.5 [2.74]	0.08 [0.03] (0.9)	0.17 [0.05] (2.0)
Fisher Urine Chemistry Control			
Level 1	61.5 [19.8]	1.1 [0.3] (1.8)	7.3 [2.3] (11.9)
Level 2	70.3 [22.7]	1.3 [0.4] (1.8)	7.3 [2.3] (10.4)
Urine Pool	43.9 [14.2]	0.6 [0.2] (1.4)	3.6 [1.2] (8.2)

e. All specific performance characteristics tests were run after normal recommended equipment quality control checks were performed. (Refer to your Dimension® Operator's Guide).

f. Specimens at each level were analyzed in triplicate for 20 days. The within-run and between-day standard deviations were calculated by the analysis of variance method.

OCS™ is a trademark of Gilford, Irvine, CA 92714.

Method Comparison

Comparative Method	Regression Statistics ^a			n
	Slope	Intercept mg/dL [mmol/L]	Correlation Coefficient	
PHOS method on the				
aca® discrete clinical analyzer				
Serum	1.00	0.0	0.976	112
Urine	1.00	-0.05 [-0.02]	0.992	122

g. Model equation for regression statistics is: [results of Dimension® system] = slope x [comparative method results] + intercept.

Specificity

HIL Interference

The PHOS method (using the 340 – 700 nm) was evaluated for interference from hemolysis, icterus and lipemia according to CLSI/NCCLS EP7-P. Bias, defined as the difference between the control sample (does not contain interferent) and the test sample (contains the interferent), is shown in the table below. Bias exceeding 10% is considered "interference".

Substance Tested	Test Concentration SI Units	PHOS Concentration mg/dL [mmol/L]	Bias % ^h
Hemoglobin (hemolysate)	50 mg/dL	1.1	<10
	[0.03 mmol/L] (monomer)	[0.36 L]	
Bilirubin (unconjugated)	5 mg/dL	1.0	<10
	[86 µmol/L]	[0.32]	
Lipemia (Intralipid®)	200 mg/dL	1.1	<10
	[2.29 mmol/L]	[0.36]	

h. Analyte results should not be corrected based on this bias.

Non-Interfering Substances

The following substances at the concentrations indicated have no significant effect (less than 10%) when added to a serum pool containing 4.0 mg/dL [1.29 mmol/L] of PHOS.

Substance	Test Concentration	SI Units
Acetaminophen	0.025 mg/dL	1.66 µmol/L
Amikacin	15 mg/dL	256 µmol/L
Ampicillin	5.3 mg/dL	152 µmol/L
Ascorbic Acid	5 mg/dL	227 µmol/L
Caffeine	6 mg/dL	308 µmol/L
Carbamazepine	3 mg/dL	127 µmol/L
Chloramphenicol	5 mg/dL	155 µmol/L
Chlordiazepoxide	1 mg/dL	33.3 µmol/L
Chlorpromazine	0.2 mg/dL	6.27 µmol/L
Cholesterol	500 mg/dL	12.9 mmol/L
Cimetidine	2 mg/dL	79.2 µmol/L
Creatinine	30 mg/dL	2652 µmol/L
Dextran 40	6000 mg/dL	1500 µmol/L
Diazepam	0.5 mg/dL	17.6 µmol/L
Erythromycin	6 mg/dL	81.6 µmol/L
Ethanol	400 mg/dL	86.8 mmol/L
Ethosuximide	25 mg/dL	1770 µmol/L
Furosemide	6 mg/dL	181 µmol/L
Gentamicin	12 mg/dL	251 µmol/L
Heparin	3 U/mL	3000 U/L
Ibuprofen	50 mg/dL	2425 µmol/L
Immunoglobulin G	5 g/dL	50 g/L
Lidocaine	1.2 mg/dL	51.2 µmol/L
Lithium	2.3 mg/dL	3.2 mmol/L
Nicotine	0.1 mg/dL	6.2 µmol/L
Penicillin G	25 U/mL	25000 U/L
Pentobarbital	8 mg/dL	354 µmol/L
Phenobarbital	10 mg/dL	421 µmol/L
Phenytoin	5 mg/dL	198 µmol/L
Primidone	4 mg/dL	183 µmol/L
Propoxyphene	0.2 mg/dL	4.91 µmol/L
Protein: Albumin	6 g/dL	60 g/L
Protein: Total	12 g/dL	120 g/L
Salicylic Acid	60 mg/dL	4.34 mmol/L
Theophylline	4 mg/dL	222 µmol/L
Urea	500 mg/dL	83.3 mmol/L
Uric Acid	20 mg/dL	1190 µmol/L
Valproic Acid	50 mg/dL	3467 µmol/L


Analytical Sensitivity: 0.0 mg/dL [0.00 mmol/L]

The analytical sensitivity represents the lowest concentration of phosphorus that can be distinguished from zero. This sensitivity is defined as the mean value (n = 20) plus two standard deviations of level 1 CHEM II calibrator.

Symbols Key: See adjacent panel.

Dimension®, Flex® and aca® are trademarks of Siemens Healthcare Diagnostics Inc.

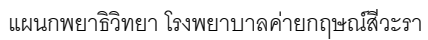©2008 Siemens Healthcare Diagnostics Inc.
All rights reserved.

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